



# **USER MANUAL**

## ***EMC-PEMF* WELLNESS MACHINE**

**Model Name: MAKINA**

**Model No. : MK2000**

## **DISCLAIMER AND CONTRAINDICATIONS**

EMC-PEMF, LLC and its distributors offer low risk, general wellness products that promote a healthy lifestyle, encourage and maintain good health, and reduce the impact or risk of some chronic diseases and conditions where a healthy lifestyle has been shown to play an important role in improving quality of life and living well. The FDA has not evaluated these statements, and EMC-PEMF products are not intended to diagnose, treat, cure or prevent any disease or medical condition. PEMF products are meant to address energetic elements that aid the body's natural healing abilities, and help to restore and maintain energetic balance. The therapeutic low intensity pulsed electromagnetic fields generated by these products use similar levels of energy as the natural vitalizing energy of the Earth. The use of a PEMF product is contraindicated in people with pacemakers or other electronic implants, cochlear implants or mechanical heart valves. PEMF products are also contraindicated during pregnancy, and in people with epilepsy or bleeding. PEMF therapy is also contraindicated for: organ transplant recipients or those on immunosuppressive therapy, and those with non-MRI safe IUD's. Patients should tell their doctors when they are thinking about using complementary therapies with conventional medical treatment. EMC-PEMF, LLC makes no warranty or representation of any kind, express or implied, regarding the suitability or fitness of any EMC-PEMF, LLC products as alternative or complementary to any medical treatments and therapies, and thus EMC-PEMF, LLC disclaims liability in connection therewith, either through the use of the purchaser or third parties. In signing this form, the purchaser acknowledges that they have been properly advised by the Terms of Use through the EMC-PEMF User Guide and Manual, and to undertake due diligence with respect to the use of the product technology. The purchaser hereby expressly acknowledges and agrees that EMC-PEMF, LLC shall not be liable for damages, or any other loss, whether direct, indirect, consequential, incidental, or special, which the purchaser may suffer, arising from any error, fault, or failure to properly use the product in accordance with the aforementioned User Guide and Manual. This information is not a substitute for professional medical advice. If you have or suspect a medical problem, please consult your health care provider.



## Table of Contents

BOX CONTENT .....	3
OPERATING INSTRUCTIONS .....	4
SET-UP PROCEDURE .....	4
USER GUIDE .....	5
CONTROL PANEL.....	6
IMPORTANT NOTES .....	8
TECHNICAL SPECIFICATIONS .....	9
TROUBLESHOOTING GUIDE.....	10
CARE AND MAINTENANCE .....	10
WARRANTY .....	11
DISCLAIMER AND CONTRAINDICATIONS.....	12
WARRANTY CARD AND NOTICE OF DISCLAIMER .....	13








A significant step towards a healthier lifestyle

Thank you for your purchase of **EMC-PEMF** Wellness  
Machine







## BOX CONTENT



Your EMC-PEMF Wellness Machine package comes with the following items.

<ul style="list-style-type: none"><li>1 pc- <i>EMC-PEMF</i> Wellness machine</li></ul>	
<ul style="list-style-type: none"><li>2 pcs- Halo Ring applicator</li></ul>	
<ul style="list-style-type: none"><li>2 pcs- 6 foot self-locking applicator cable</li></ul>	
<ul style="list-style-type: none"><li>1 pc- 6 foot two-prong AC cord</li></ul>	
<ul style="list-style-type: none"><li>1 pc- Waterproof accessory hard-case</li></ul>	
<ul style="list-style-type: none"><li>1 pc – <i>EMC-PEMF</i> authenticator neodymium disc</li></ul>	
<ul style="list-style-type: none"><li>1 pc- 2-ft x 6-ft foldable <i>EMC-PEMF</i> Gaia Mat (optional accessory upgrade)</li></ul>	

## OPERATING INSTRUCTIONS

### SET-UP PROCEDURE

<p>A. Insert AC cord to the AC inlet, then plug in the machine to AC supply outlet</p> <p><i>EMC-PEMF</i> wellness machine will operate with an AC supply of 100 VAC – 240 VAC 50/60Hz</p>	
<p>B. Gently connect applicator cable to the HALO RING</p>	
<p>C. Gently insert applicator cable/s into either port P1 or P2, or both ports of the machine</p>	
<p>D. Press the power button to turn-on</p>	

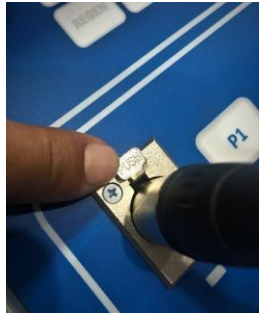
<p>E. <b>Caution:</b> DO NOT forcefully pull out applicator cable when unplugging from the port to avoid damage. Push release button, then pull out applicator cable</p>	
<p>F. <b>Caution:</b> DO NOT forcefully pull out applicator cable when disconnecting the applicator ring to avoid damage. Push release button, then pull out applicator cable</p>	

## USER GUIDE

After set-up, the machine is now ready to operate. Upon pressing the power button to turn-on, the panel light indicators are default to **MODE- RELAX** setting.



## CONTROL PANEL

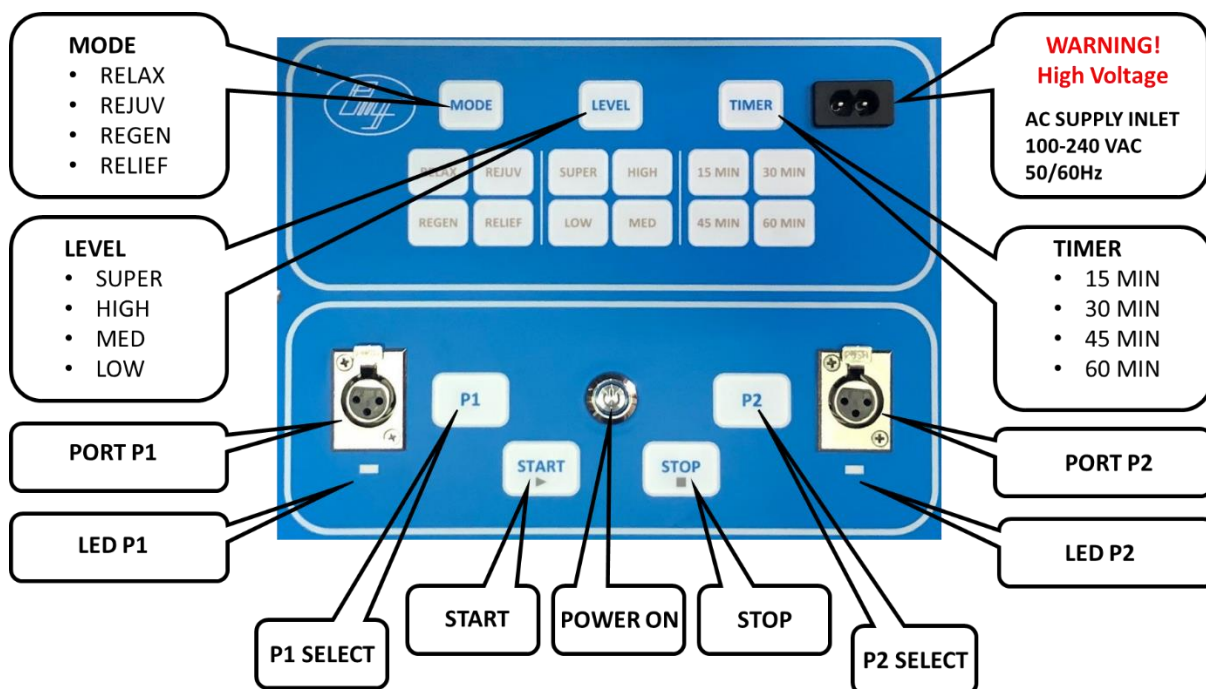

- 1- Press **MODE** button to select one of the preset programs. Continue pressing until desired program lights up. There are four **MODE** preset programs to choose from: **RELAX, REJUV, REGEN, RELIEF**



To override default setting of **LEVEL** or **TIMER** continue to steps two (2) and three (3). Otherwise, continue to step four (4).


### MODE Preset Program Selection

MODE (Program)	LEVEL (PEMF)	TIMER (Minute)
RELAX	LOW	60 MIN
REJUV	MED	15 MIN
REGEN	MED	30 MIN
RELIEF	HIGH	30 MIN


- 2- Press **LEVEL** button to select intensity level of the magnetic field. Continue pressing until desired program lights up. There are four intensity levels to select from: **LOW**, **MED**, **HIGH**, **SUPER** (i.e. maximum)



- 3- Press **TIMER** button to select the length of time to operate the machine. Continue pressing until desired program lights up. There are four timer settings to choose from: **15 MIN**, **30 MIN**, **45 MIN**, **60 MIN** (MIN – minute)

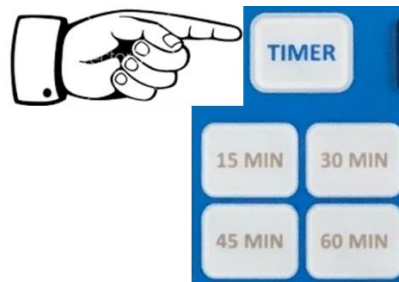

- 4- Press either port **P1** or port **P2** button, or both to operate two ring applicators. When pressed the indicator will light up.

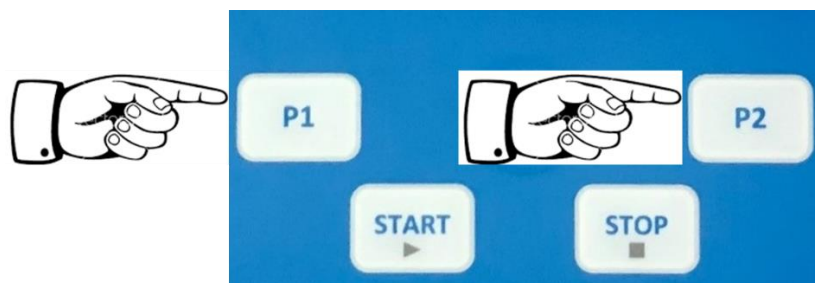



- 5- Press **START** button to run the program selected. P1 and/or P2 light indicator will start blinking, which indicates the program is running. **START** button will not respond if neither port **P1** nor **P2** is on. Make sure to press either port **P1** or **P2** button or both before pressing **START** to run program.



- 6- Press **STOP** button more than one second to end program anytime. Otherwise, the machine will run program until the selected timer ends. A long consecutive continuous beeping sound signifies the end of the program and **P1** and/or **P2** light will turn off.



## IMPORTANT NOTES

- a- After pressing **START**, changing **MODE**, **TIMER**, and **LEVEL**, **P1** or **P2** setting will not be possible. To make changes press the **STOP** button.
- b- Pressing **STOP** button will always reset **TIMER**.
- c- **START** button will not respond if neither port **P1** nor **P2** is on. Make sure to press either port **P1** or **P2** button or both before pressing **START** to run program.

## TECHNICAL SPECIFICATIONS

Model Name	MAKINA	
Model Number	MK2000	
Features		
<ul style="list-style-type: none"><li>Selectable Mode</li><li>Selectable Intensity Level</li><li>Selectable Timer</li><li>Applicator Ports</li><li>Modulation Pattern</li><li>Intensity Range</li><li>Magnetic Field Pattern</li><li>Light Indicator</li></ul>	<ul style="list-style-type: none"><li>➤ 4- (Relax, Rejuv, Regen, Relief)</li><li>➤ 4- (Super (Max), High, Med, Low)</li><li>➤ 4- (15, 30, 45, 60 minutes)</li><li>➤ 2- (P1, P2)</li><li>➤ Triangular PWM</li><li>➤ 0 – 10 milliTesla</li><li>➤ Pulsed Quasi-Sinusoid</li><li>➤ Blue LED</li></ul>	
	Parameter	Unit
Electrical Specifications		
AC Input Voltage	100 - 240	Vac
AC Input Current	1.66	Amp
Frequency	50/60	Hz
Peak Power	200	Watt
Mechanical Specifications		
Operating Temperature	0° to 35°	Celsius
Storage Temperature	-15° to 60°	Celsius
Case Dimension	12.5× 10.1 × 6.0”	(L x W x H): in
Approximate Weight	4.5	kg
Environment	Waterproof (IPX7), dustproof (IP6X) & MIL-SPEC	



## TROUBLESHOOTING GUIDE

Issue	Check if...	Then
Control panel does not respond and/or LED indicator lights are off	AC cord is plugged to the AC inlet of the machine and AC cord plugged to AC supply	Press power button
Applicator ring is not pulsing magnetic field, but LED indicator is flashing	Applicator cable is properly inserted to the port (P1 or P2) and applicator ring	Hold neodymium disc in the middle of the ring to feel pulsing energy
Applicator ring is not pulsing magnetic field, but port LED indicator is OFF	Port P1 or P2 is pressed	Press Start button. Hold neodymium disc in the middle of the ring to feel pulsing energy
Applicator ring is not pulsing magnetic field, but port LED indicator is ON steady	Start button is pressed	Press Start button. Hold neodymium disc in the middle of the ring to feel pulsing energy

## CARE AND MAINTENANCE

ITEM	SUGGESTIONS
MAKINA DEVICE	Keep out of reach of children; Keep dry; Do not use in or around water (i.e. bathtub, sink, pool); Do not tamper; Wipe clean after use
CASE	Do not drop; Do not spill water when it is open; Store within recommended temperature range; Wipe off dust
HALO RING	Do not submerge in water; Do not bend ring; Do not forcefully pull connector, use release latch to disconnect cable; Do not peel off protective cover; Wipe clean after use
CONTROL PANEL	Keep dry; Avoid water/liquid spill; Use release latch to disconnect applicator cable; Gently pull AC cord out of the AC inlet; Wipe clean after use
APPLICATOR CABLE	Keep dry; Do not submerge in water; Do not forcefully pull connector, use release latch to disconnect cable; Wipe clean after use

For any questions or clarifications, please message [customercare@emcpemf.com](mailto:customercare@emcpemf.com), and we'll get back to you within seven (7) business days.



## WARRANTY

**Your complete satisfaction is important to us.** EMC-PEMF, LLC warrants each new EMC-PEMF Wellness machine, EMC-PEMF Halo Ring, EMC-PEMF Venus Ring, EMC-PEMF Gaia-Mat, EMC-PEMF ACCORDION, and EMC-PEMF D'ROD to be free from defects in materials and workmanship for a period of THREE (3) YEARS after the sale of the EMC-PEMF product to the original purchaser. The obligation of EMC-PEMF under this warranty is expressly limited solely and exclusively to the repair or replacement of the unit or any parts thereof, which to EMC-PEMF satisfaction shall have become defective during the warranty period, and which shall have been returned to EMC-PEMF within five (5) business days after the discovery of the defect by the original purchaser. This warranty does not extend to any liability for medical expenses, or for any other defect or malfunction of the EMC-PEMF product, except as herein provided, whether such damage claim shall be based on contract, tort, breach of warranty, or otherwise. This warranty shall not apply to any EMC-PEMF product which has been repaired, tampered with or altered by someone other than a duly authorized EMC-PEMF representative, or which has been subjected to negligence, accident, mishandling, normal wear and tear or which has not been used in accordance with the enclosed instructions or for the stated purposes. This warranty is expressly limited solely to the original purchaser (user) and does not extend to any transferee, assignee, or subsequent purchaser or user of the EMC-PEMF product. This warranty is the sole and exclusive warranty made or which may be deemed to have been made by EMC-PEMF and is expressly in lieu of any and all other warranties, either expressed or implied, including any warranties of merchantability or fitness for a particular purpose. No person or entity has any authority to bind EMC-PEMF to any warranty, guarantee or representation except as specifically set forth herein.

PRODUCT SERIAL NUMBER	
ORDER NUMBER	
CLIENT NAME & SIGNATURE	
DATE OF PURCHASE	



## DISCLAIMER AND CONTRAINDICATIONS

EMC-PEMF, LLC and its distributors offer low risk, general wellness products that promote a healthy lifestyle, encourage and maintain good health, and reduce the impact or risk of some chronic diseases and conditions where a healthy lifestyle has been shown to play an important role in improving quality of life and living well. The FDA has not evaluated these statements, and EMC-PEMF products are not intended to diagnose, treat, cure or prevent any disease or medical condition. PEMF products are meant to address energetic elements that aid the body's natural healing abilities, and help to restore and maintain energetic balance. The therapeutic low intensity pulsed electromagnetic fields generated by these products use similar levels of energy as the natural vitalizing energy of the Earth. The use of a PEMF product is contraindicated in people with pacemakers or other electronic implants, cochlear implants or mechanical heart valves. PEMF products are also contraindicated during pregnancy, and in people with epilepsy or bleeding. PEMF therapy is also contraindicated for: organ transplant recipients or those on immunosuppressive therapy, and those with non-MRI safe IUD's. Patients should tell their doctors when they are thinking about using complementary therapies with conventional medical treatment. EMC-PEMF, LLC makes no warranty or representation of any kind, express or implied, regarding the suitability or fitness of any EMC-PEMF, LLC products as alternative or complementary to any medical treatments and therapies, and thus EMC-PEMF, LLC disclaims liability in connection therewith, either through the use of the purchaser or third parties. In signing this form, the purchaser acknowledges that they have been properly advised by the Terms of Use through the EMC-PEMF User Guide and Manual, and to undertake due diligence with respect to the use of the product technology. The purchaser hereby expressly acknowledges and agrees that EMC-PEMF, LLC shall not be liable for damages, or any other loss, whether direct, indirect, consequential, incidental, or special, which the purchaser may suffer, arising from any error, fault, or failure to properly use the product in accordance with the aforementioned User Guide and Manual. This information is not a substitute for professional medical advice. If you have or suspect a medical problem, please consult your health care provider.



===== CUT HERE =====

## WARRANTY CARD AND NOTICE OF DISCLAIMER

The USER MANUAL, OPERATING INSTRUCTIONS, WARRANTY, DISCLAIMER, and CONTRAINDICATIONS form part of the “CLIENT–EMC-PEMF CONTRACT OF WARRANTY AND WAIVER OF LIABILITIES FOR IMPROPER USE AND HANDLING.” Please read them with care. For purposes of PURCHASER CARE, QUALITY ASSURANCE and WARRANTY, as well as AGREEMENT to the PROPER TERMS OF USE and PRODUCT DISCLAIMERS therein, kindly affix your name, signature and date on the space provided below. This serves as your ACKNOWLEDGEMENT and CONSENT regarding limited product and medical liability, return and repair in case of manufacturer’s defect, terms of use, and the full warranty period of THREE (3) YEARS, among others, as outlined in the USER MANUAL, OPERATING INSTRUCTIONS, WARRANTY, DISCLAIMER, and CONTRAINDICATIONS.

PRODUCT SERIAL NUMBER	
ORDER NUMBER	
CLIENT NAME & SIGNATURE	
DATE OF PURCHASE	

Please send completed form to:

**EMC-PEMF, LLC**  
**497 Suisse Drive**  
**San Jose, CA 95123, Silicon Valley**  
**United States**

Alternatively, please email scanned signed copy to:

[warranty@emcpemf.com](mailto:warranty@emcpemf.com)

You will receive a notification of receipt. This signed document may be used as LEGAL BASIS in case of loss of the WARRANTY, DISCLAIMER, and CONTRAINDICATIONS, upon claim of PRODUCT LIABILITY and WARRANTY.

